



# Morgan™ HLA SSP ABC Typing Kit

## Instructions for Use

**REF** 33210

**REF** 33211

**IVD**

**CE** 0197

Sep. 2014  
Rev. 4.0

## Table of Contents

<b>1. INTENDED USE .....</b>	<b>2</b>
<b>2. SUMMARY AND EXPLANATION .....</b>	<b>2</b>
<b>3. PRINCIPLE(S) .....</b>	<b>3</b>
<b>4. REAGENTS .....</b>	<b>4</b>
4.1 Contents of the Morgan™ HLA SSP ABC Typing Kit	
4.2 Warning or Caution	
4.3 Storage	
4.4 Instability Indications	
<b>5. INSTRUMENT REQUIREMENTS .....</b>	<b>6</b>
5.1 Programming the Thermal Cycler	
5.2 Gel Electrophoresis	
<b>6. SPECIMEN COLLECTION AND PREPARATION .....</b>	<b>6</b>
6.1 Blood Sample Collection	
6.2 DNA Isolation	
6.3 DNA Quantity	
6.4 DNA Quality	
<b>7. DIRECTIONS FOR USE .....</b>	<b>7</b>
7.1 Materials Provided	
7.2 Materials Required, But Not Provided	
7.3 Sample Preparation	
7.4 Reagent/Equipment Preparation	
7.5 Performing PCR	
7.6 Gel Electrophoresis	
<b>8. RESULTS .....</b>	<b>10</b>
8.1 Data Analysis	
8.2 Gel Interpretation	
8.3 Limitation of Procedure	
8.4 Quality Control	
<b>9. PERFORMANCE CHARACTERISTICS .....</b>	<b>12</b>
9.1 Specific Performance	
9.2 Quality Control of Kit Manufacturing	
9.3 Reproducibility	
<b>10. TROUBLESHOOTING .....</b>	<b>13</b>
<b>11. BIBLIOGRAPHY .....</b>	<b>14</b>
<b>12. TRADEMARKS USED IN THIS DOCUMENT .....</b>	<b>14</b>
<b>13. PATENTS USED IN THIS DOCUMENT .....</b>	<b>14</b>

## 1. INTENDED USE

This kit provides the HLA-A, HLA-B and HLA-C typing results with a low to medium resolution by a DNA-based method.

## 2. SUMMARY AND EXPLANATION

The HLA molecules of the human major histocompatibility (MHC) antigens are key elements in restricting the specificity of T-cell mediated immune responses.<sup>1</sup> As HLA compatibility between donor and recipient is a major factor for determining the outcome of organ and bone marrow transplantations, accurate determination of HLA genotypes is necessary prior to transplantation. Traditionally, HLA typing has been performed by serological methods.<sup>2</sup> However, limitations of the traditional serological methodology become apparent as we become more knowledgeable of the HLA genes.

The HLA genes, located on the short arm of human chromosome 6, are highly polymorphic. As more DNA sequences of HLA alleles are revealed, it was clear that most HLA antigens differ from each other by only a few or one single amino acid. Due to the limitation of the traditional serological methodology differentiate these subtle antigenic differences, the use of DNA-based methods were introduced to provide a higher typing resolution. TBG Biotechnology would like to introduce the Morgan<sup>TM</sup> HLA SSP ABC Typing Kit as a DNA-based method for the determination of HLA alleles using sequence specific primers (SSP).<sup>3</sup>

<p><b>WARNING:</b> The results obtained from the Morgan<sup>TM</sup> HLA SSP ABC Typing Kit should not be used as the sole basis for making a clinical decision.</p>
--

### **3. PRINCIPLE(S)**

The Morgan™ HLA SSP ABC Typing Kit is for determining HLA alleles using PCR techniques with sequence specific primers (PCR-SSP). Allele sequence-specific primer pairs are designed to selectively amplify target sequences that are specific to a single allele or group of alleles. This PCR-SSP method is based on the principle that only primers with completely matched sequences to the target sequences result in amplified products under controlled PCR conditions. The presence of amplified DNA fragment is a positive indication of the existence of allele specific sequence within the genomic DNA. On the other hand, mismatched primers do not generate amplicons.

In addition to sequence specific primers, an internal control primer pair, which amplifies a conserved region of the house keeping gene, cystic fibrosis gene, is included in every PCR reaction mix, and the PCR product of the internal control primer pair serves as an indication of the integrity of PCR reaction. When the DNA fragment amplified by sequence specific primer is present, the product of the internal control primer pair is sometimes weak or absent due to the differences in amplifying efficiency between the specific primer pair and the internal control primer pair.

The PCR reaction products are examined by agarose gel electrophoresis that separates the DNA fragments by size. The amplified fragments are visualized as bands through ethidium bromide staining followed by ultraviolet light exposure.

The Morgan™ HLA SSP ABC Typing Kit supplies a panel of primer mixes that is designed to provide positive identification of HLA alleles. A worksheet is provided to assist the determination of HLA types by analyzing the resulting PCR reactions.


## 4. REAGENTS

### 4.1 Contents of the Morgan™ HLA SSP ABC Typing Kit

The Morgan™ HLA SSP ABC Typing Kit provides sufficient reagents for 10 tests or 20 tests.

#### 4.1.1 Typing Trays

Each tray is designed for 1 test, consisting of PCR tubes that contain dried primers. A black dot marks well one, which is also the contamination control well, for each test. Tray layout with primer mix number is shown below.

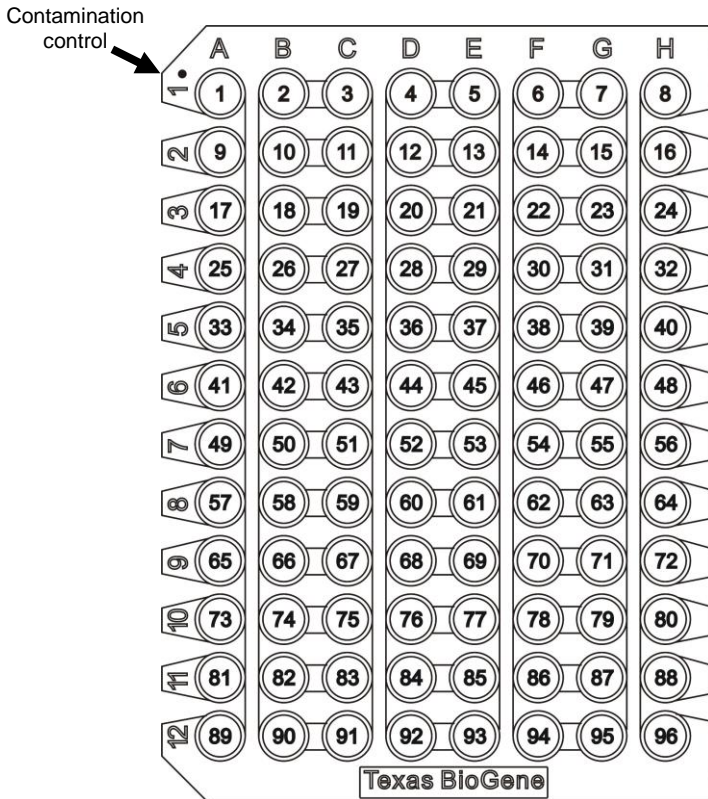

Contamination Control: 1

HLA-A: 2-24

HLA-B: 25-72

HLA-C: 73-96

#### **4.1.2 Master Buffer 1020µl**

The buffer is optimized for PCR reactions, and contains nucleotides, salts, glycerol and cresol red (final concentration: 0.15 mM dNTP, 1.5 mM MgCl<sub>2</sub>, 50 mM KCl, 15.75 mM Tris-HCl, 0.5 % glycerol). Each pre-aliquoted tube is for 1 test.

#### **4.1.3 PCR Plate Cover Membrane**

The cover membrane is used to prevent evaporation loss during the PCR process.

#### **4.1.4 Gel Interpretation Worksheets**

The worksheet reveals the primer sets that are reactive to specific HLA alleles.

### **4.2 Warning or Caution**

For In Vitro Diagnostic Use

**4.2.1** Do not use reagents past the expiration date printed on the label.

**4.2.2** The test must be performed by well-trained and authorized laboratory staff.

**4.2.3** Reagents should be handled in accordance with good laboratory practice using appropriate precautions.

**4.2.4** Pipettes used for post-PCR manipulations should not be used for pre-PCR manipulations.

**4.2.5** Biohazard Warning: the ethidium bromide used for staining of DNA is a potential carcinogen. Always wear gloves when handling stained gels.

**4.2.6** Biohazard Warning: all blood products should be treated as potentially infectious.

**4.2.7** Caution: wear UV-blocking eye protection and do not view UV light source directly when viewing or photographing gels.

See Material Safety Data Sheets for detailed information.

### **4.3 Storage**

The Typing Trays are sealed in pouches and should be stored at 2-8°C. The Master Buffer should be stored at -20°C. Any unused Typing Tray should be returned to their original pouch and sealed using the zipper. Once the pouch is opened, the Typing Tray should be used within 4 weeks.

### **4.4 Instability Indications**

Do not use Typing Trays with cracks in the tubes.

## 5. INSTRUMENT REQUIREMENTS

### 5.1 Programming the Thermal Cycler

The Morgan™ HLA SSP ABC Typing Kit is standardized for the following program using Applied Biosystems GeneAmp® 9600 and 9700. Thermal cyclers other than the recommended model and/or brand have to be user-validated. The program runs in approximately 1 hr 25 min.

**Morgan HLA Program**

Segment	Cycle Number	Temperature	Time
1	1	96°C	2.5 min
2	10	96°C	15 sec
		65°C	60 sec
3	22	95°C	15 sec
		62°C	50 sec
		72°C	30 sec
4	1	4°C	until removed

Please program the thermal cycler according to manufacture's instruction.

### 5.2 Gel Electrophoresis

Refer to Section 7.6 Gel Electrophoresis.

## 6. SPECIMEN COLLECTION AND PREPARATION

### 6.1 Blood Sample Collection

The peripheral blood can be collected using blood collection tubes with sodium citrate and EDTA as anticoagulants. Contamination of the DNA by heparin can result in interference with the PCR reaction. For this reason, heparinized blood should not be used as a starting material for DNA isolation. Please refer to manufacturer's instruction for blood sample storage and stability.

### 6.2 DNA Isolation

Genomic DNA may be obtained from all nucleated cells. The source of cell includes, but is not limited to, whole blood, buffy coat and cultured cells. DNA isolation can be performed by any validated protocol that produces DNA with adequate quality and quantity of DNA followed by repeated user validation.

### 6.3 DNA Quantity

The DNA sample should be re-suspended in sterile distilled water or appropriate buffer solution at a concentration of 10-40 ng/μl. DNA should

not be re-suspended in solutions containing chelating agents, such as EDTA, above 0.5 mM in concentration.

#### **6.4 DNA Quality**

The DNA sample should have an A260/A280 ratio between 1.65 and 1.8.

DNA samples may be used immediately after isolation or stored at  $-20^{\circ}\text{C}$  or below for extended periods of time (over 1 year) with no adverse effects on the HLA typing results.

## **7. DIRECTIONS FOR USE**

### **7.1 Materials Provided**

Refer to Section 4. REAGENTS

### **7.2 Materials Required, But Not Provided**

#### **7.2.1 PCR**

- 7.2.1.1 Taq DNA polymerase (5 U/ $\mu\text{l}$ ) (e.g. Fermentas, Cat. #EP0401 or EP0402)
- 7.2.1.2 Distilled water or DNA diluent
- 7.2.1.3 Adjustable pipettes
- 7.2.1.4 Disposable pipette tips
- 7.2.1.5 Vortex mixer
- 7.2.1.6 Microcentrifuge
- 7.2.1.7 Table top centrifuge for 96 well PCR plates
- 7.2.1.8 Thermal Cycler with heated cover (e.g. Applied Biosystems GeneAmp 9600 and 9700)
- 7.2.1.9 Pressure pad

#### **7.2.2 Gel Electrophoresis**

- 7.2.2.1 Molecular biology grade agarose
- 7.2.2.2 0.5x TAE or TBE buffer
- 7.2.2.3 ethidium bromide solution (10 mg/ml)
- 7.2.2.4 Hot plate or microwave oven for heating agarose solutions
- 7.2.2.5 Gel tray and electrophoresis chamber Power supply
- 7.2.2.6 UV transilluminator
- 7.2.2.7 Photographic or image documentation system.

### **7.3 Sample Preparation**

Refer to Section 6. SPECIMEN COLLECTION AND PREPARATION.

### **7.4 Reagent/Equipment Preparation**

Refer to Section 5. INSTRUMENT REQUIREMENTS.



## 7.5 Performing PCR

- 7.5.1 Remove the tube of pre-aliquoted Master Buffer from freezer and thaw at room temperature (20-25°C). Each tube is for one test.
- 7.5.2 Once the Master Buffer is thawed, vortex to mix well. Pulse-spin the buffer tube in a microcentrifuge to bring all liquid down to the bottom.
- 7.5.3 Remove Typing Trays from fridge. Take out the trays from the pouch. Return any unused tray to pouch and seal with the pouch zipper for future use.
- 7.5.4 Remove Taq DNA polymerase (5 U/μl) from freezer and keep on ice.
- 7.5.5 Add 2 μl of distilled water or DNA diluent to the first well marked by a black dot. The first well is also the contamination control well.
- 7.5.6 Add 6.5 μl of Taq polymerase (5U/μl) to the Master Buffer tube. Cap tube and vortex to mix well. Pulse-spin in a microcentrifuge to bring all liquid down to the bottom.
- 7.5.7 Dispense 10 μl of the Master Buffer/Taq polymerase mixture to the contamination control well.
- 7.5.8 To make the master mix, add 202 μl of DNA sample (10-40 ng/μl) to the remaining Master Buffer/Taq polymerase mixture. Cap tube and vortex to mix well. Pulse-spin in a microcentrifuge to bring all liquid down to the bottom.
- 7.5.9 Dispense 12 μl of the master mix into wells 2 through 96 of the Typing Tray. Usage of a repeater pipette is recommended. Care should be taken to avoid cross-contamination between wells by applying the master mix to the walls of the tube. Do not allow the pipette tips to contact the dried primer pellet at the bottom of each well.
- 7.5.10 Cover the Typing Tray with PCR Plate Cover Membrane. Make sure all tubes are properly covered by the membrane to prevent evaporative loss during the PCR process. Make sure the dispensed solution comes in contact with the dried primers on bottom of each tube by tapping or centrifugation.
- 7.5.11 Place a pressure pad on top of the PCR Plate Cover Membrane before closing the thermal cycler lid.
- 7.5.12 Close the heated thermal cycler lid and start the program shown in Section 5.1 of this manual.
- 7.5.13 After the PCR process is completed, you may proceed to gel electrophoresis immediately or store the tray at -20 °C and continue to gel electrophoresis at a later time.

## **7.6 Gel Electrophoresis**

### **7.6.1 Gel Preparation**

- 7.6.1.1 Prepare sufficient 2% agarose gel to resolve 96 PCR reactions according to the original gel electrophoresis system manufacture's instructions.
- 7.6.1.2 Make appropriate adjustments to ensure that the gel contains a final ethidium bromide concentration of 0.5 µg/ml.
- 7.6.1.3 Make sure that each well of the gel creates a minimum volume of 15 µl.
- 7.6.1.4 Make sure that the running buffer used for electrophoresis is consistent with the buffer used for gel preparation.

### **7.6.2 Performing Gel Electrophoresis**












- 7.6.2.1 Transfer the entire volume of each PCR reaction mix (12 µl) into the wells in the proper sequence. A black dot marks the position one for each test and the order of samples is from left to right, top to bottom.
- 7.6.2.2 The electrophoresis takes place at 10V/cm until the cresol red has migrated about 0.7-1.0 cm into the gel.
- 7.6.2.3 After the electrophoresis is completed, place the agarose gel on an UV transilluminator. Photograph for documentation and HLA type interpretation.

## 8. Results

### 8.1 Data Analysis

- 8.1.1** The DNA fragments amplified by allele- or group- specific primer pairs are designed to have a size in the range between 70-275 bp. Refer to the worksheet for the exact size of the specifically amplified DNA fragment in each reaction tube. The internal control fragment is 600 bp in size. Therefore, the internal control band is the slower migrating band compared to the allele- or group- specific band. The presence of the allele- or group- specific band is scored as positive for the particular reaction. The absence of this allele- or group- specific band is scored as negative.
- 8.1.2** The internal control band may be weak or absent when the allele- or group- specific band is present. As the presence of the specific band is also an indication to the success of the PCR reaction. This is not a limitation to the test.
- 8.1.3** The internal control band should always be visible in negative lanes. Absence of both the internal control band and allele- or group- specific band indicates failed reaction and can void the test result.
- 8.1.4** The presence of a band in the contamination control well voids the test result.
- 8.1.5** Document the positive reactions for the test and determine the HLA types by using the worksheet provided.

### 8.2 Gel Interpretation

	Positive Reaction	Positive Reaction	Negative Reaction	Non-amplification
Well				
Control band	None	- - - - -		None
Specific band			None	None
Primer dimer				

### **8.3 Limitation of Procedure**

- 8.3.1** Performance of the test can only be guaranteed if the enclosed instructions are strictly followed.
- 8.3.2** All instruments, including the thermal cycler and pipetting devices, must be calibrated according the manufacturer's instructions.
- 8.3.3** The Morgan™ HLA SSP ABC Typing Kit is standardized for the Taq polymerase from Fermentas. The use of other brands of Taq polymerase may result in non-specific false positive PCR reactions.
- 8.3.4** The quality and quantity of the DNA samples have to be within the range as specified within instruction manual.
- 8.3.5** The Morgan™ HLA SSP ABC Typing Kit can not resolve all allele combinations. When an ambiguous result obtained from the Morgan™ HLA SSP ABC Typing Kit, please refer to section 10, TROUBLESHOOTING, to identify possible causes and solutions. If the ambiguous results are still not resolved, please use a typing kit with the same or higher level of resolution from other suppliers to resolve the ambiguity. You can also contact TBG Biotechnology Corp. for assistance.
- 8.3.6** The results obtained from the Morgan™ HLA SSP ABC Typing Kit should not be used as the sole basis for making a clinical decision.

### **8.4 Quality Control**

Each manufactured lot is checked against a panel of standard DNA samples. The alleles carried by these DNA samples react with the corresponding primer mixes in the kit. The report is available upon request.

## 9. PERFORMANCE CHARACTERISTICS

### 9.1 Specific Performance

The Morgan™ HLA SSP ABC Typing Kit was compared with two SSP and one SSO typing system by testing 252 random whole blood specimens at six geographically distinct locations. There was 100% concordance (252/252) for the HLA-A alleles and 99% concordance (250/252) for the HLA-B, and 100% concordance (252/252) for the HLA-C alleles between the Morgan™ HLA SSP ABC Typing Kit results and results obtained using the other SSP and SSO products, as shown in the following table.

Allele	Number of Concordant Tests	Number of Tests	Concordance	95 % Confidence Limits	
				Lower	Upper
HLA-A	252	252	100%	98.8%	100%
HLA-B	250	252	99.2% <sup>a</sup>	97.5%	100%
HLA-C	252	252	100% <sup>b</sup>	98.85%	100%

- a) One of the two non-concordant results for the B alleles is a failed test which might be due to DNA specimen problem as suggested by the off-site study director. The other one is an ambiguous result. One of the options contained within the ambiguous results matched the unambiguous results of the predicated product.

### 9.2 Quality Control of Kit Manufacturing

Each manufactured lot is checked against a panel of standard DNA samples. The alleles carried by these DNA samples react with most of the primer mixes in the kit. The report is available upon request.

### 9.3 Reproducibility

Ten standard DNA samples from International Histocompatibility Working Group (IHWG) were typed using the Morgan™ HLA SSP ABC Typing Kit for 5 times for the within-lot and inter-lot reproducibility studies. The data show 100% reproducibility for HLA-A, HLA-B and HLA-C alleles.

## 10. TROUBLESHOOTING

<b>Problem</b>	<b>Potential Cause</b>	<b>Solution</b>
Weak bands or absence of bands	Poor DNA quality	Repeat DNA purification and make sure the A260/A280 ratio is between 1.65-1.80.
	Insufficient amount of DNA	Make sure the DNA concentration is within the suggested range of 10~40ng/μl.
	Insufficient Taq polymerase activity	Make sure there is sufficient amount of Taq polymerase added and the Taq polymerase has not been degraded.
	Improperly calibrated thermal cycler	Calibrate the thermal cycler.
	Unevenly mixed PCR master mix	Make sure the master mix is thoroughly mixed before use.
	Insufficient ethidium bromide staining	Re-stain the agarose gel in fresh 0.5XTBE buffer with 0.5 μg/ml ethidium bromide.
Presence of false positive band(s)	Overly high DNA concentration	Make sure the DNA concentration is within the suggested range of 10~40ng/μl.
	Excessive Taq polymerase activity	Make sure the amount of Taq polymerase added is as instructed on page 8.
	Improperly calibrated thermal cycler	Calibrate the thermal cycler.
Band(s) in the contamination control well	Reagent contamination	Use new reagents for PCR preparation.

## **11. BIBLIOGRAPHY**

1. Cruse JM and Lewis RE, 1998, Atlas of Immunology. CRC Press: 77-97
2. Terasaki, P.I., Bernoco, F., Park, M.S., Ozturk, G., and Iwaki, Y. Microdroplet testing for HLA-A, -B, -C and -D antigens. American Journal of Clinical Pathology 69:103-120, 1978.
3. Newton, C.R., Graham, A., Heptinstall, E., Powell, S.J., Summers, C., Kalsheker, N., Smith, J.C., and Markham, A.F. Analysis of any point mutation in DNA. The amplification refractory mutation system (ARMS). Nucleic Acids Research 17:2503-2516, 1989.

## **12. TRADEMARKS USED IN THIS DOCUMENT**

Morgan™  
TBG Biotechnology Corp.  
Applied Biosystems  
GeneAmp®

## **13. PATENTS USED IN THIS DOCUMENT**

This product is optimized for use in the Polymerase Chain Reaction (“PCR”) Process which is covered by patents owned by Roche Molecular Systems, Inc. and F. Hoffmann-La Roche Ltd (“Roche”). No license under these patents to use the PCR Process is conveyed expressly or by implication to the purchaser by the purchase of this product. Further information on purchasing licenses to practice PCR may be obtained by contacting, in the United States, the Director of Licensing at Roche Molecular Systems, Inc. 1145 Atlantic Avenue, Alameda, California 94501, and outside the United States, the PCR Licensing Manager, F. Hoffmann-La Roche Ltd, Grenzacherstr. 124, CH-4070 Basel, Switzerland.

### **European Representative:**

mdi Europa GmbH, Langenhagener Str. 71, 30855 Hannover-Langenhagen, Germany, +49 511 39 08 95 30



TBG Biotechnology Corp., 13F., No.237, Sec. 1, Datong Rd.,  
Xizhi Dist., New Taipei City 221, Taiwan (R.O.C.)